

Ramic

- Intramedullary Nail Systems
- Humeral Nail



Features




Straight nail with central insertion point, to avoid potential insertion through the fracture site, and to reduce the harm on supraspinatus tendon.

Distal locking with 3.5mm screw:
Short nail, interlocking design to provide more stability.
Long nail, special screw angle to avoid the radial nerve and provide secured locking position.

Three size end cap are available for options.

Proximal locking with 4.0mm screw: left and right design to provide stable fixation on Humeral Head . Ascending screw hole to avoid the harm on allxiary nerve and posterior humeral circumflex artery.

Blunt screw tip to reduce harm on soft tissues. Low profile head to minimize the screw exposure .

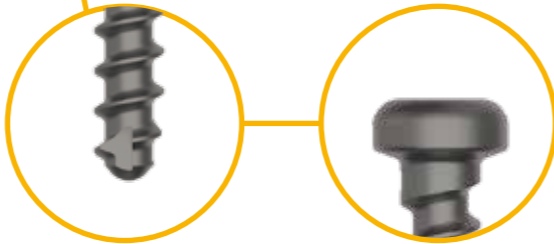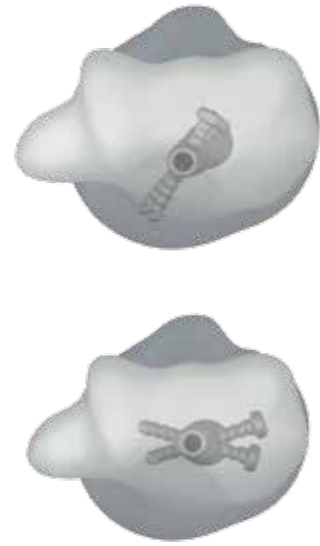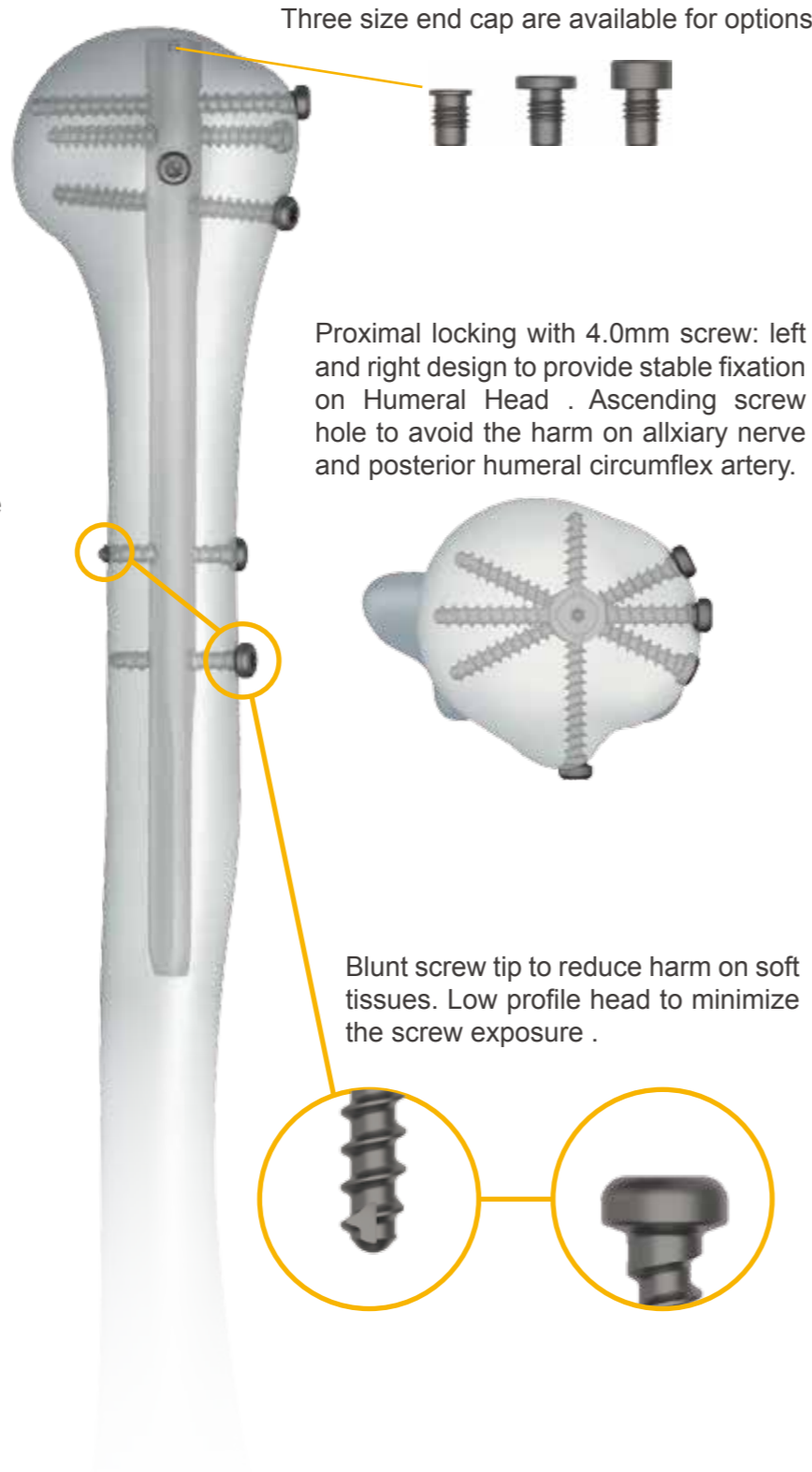

Humeral Nail-Left

Product Code	Diameter	Length
55407016	Φ7	160
55407020	Φ7	200
55407022	Φ7	220
55407024	Φ7	240
55407026	Φ7	260
55407028	Φ7	280
55407030	Φ7	300
55408016	Φ8	160
55408020	Φ8	200
55408022	Φ8	220
55408024	Φ8	240
55408026	Φ8	260
55408028	Φ8	280
55408030	Φ8	300
55409016	Φ9	160

Humeral Nail-Right

Product Code	Diameter	Length
55417016	Φ7	160
55417020	Φ7	200
55417022	Φ7	220
55417024	Φ7	240
55417026	Φ7	260
55417028	Φ7	280
55417030	Φ7	300
55418016	Φ8	160
55418020	Φ8	200
55418022	Φ8	220
55418024	Φ8	240
55418026	Φ8	260
55418028	Φ8	280
55418030	Φ8	300
55419016	Φ9	160



Humeral Locking Screw - Φ3.5

Product Code	Length
55423520	20
55423522	22
55423524	24
55423526	26
55423528	28
55423530	30
55423532	32
55423534	34
55423536	36
55423538	38
55423540	40
55423542	42
55423544	44
55423546	46
55423548	48
55423550	50
55423552	52
55423554	54
55423556	56
55423558	58
55423560	60



Humeral Locking Screw - Φ4.0

Product Code	Length
55424022	22
55424024	24
55424026	26
55424028	28
55424030	30
55424032	32
55424034	34
55424036	36
55424038	38
55424040	40
55424042	42
55424044	44
55424046	46
55424048	48
55424050	50
55424052	52
55424054	54
55424056	56
55424058	58
55424060	60



Indications

Proximal humeral fractures Humeral shaft fractures

Humeral End Caps

Product Code	Extension
55430000	0
55430002	2
55430004	4

